

Competitor II Semi Electric Height Adjustable Bed

HCPCS:



FIGURE A



FIGURE B



FIGURE C



FIGURE D

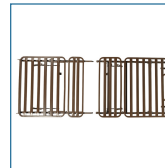


FIGURE E

- #15571 Hand Crank adjustable bed height features a NEW positive fit hexagonal crank connection making it easy to raise and lower bed height for patient comfort, ease of entry/egress and patient care accessibility
- Large glow in the dark hand pendant has easy-to-use controls for adjusting head and foot deck sections to comfortable positions built in battery backup and power indicator light. NEW Lock out feature can be used to prevent accidental actuation of head or foot adjustment
- In the event of a power failure, a 9 volt battery (not included) located in the hand pendant can lower the head and foot sections to a horizontal position. Two additional Batteries can be installed into the motor housing to increase adjustment speed and extend battery life. No crank needed
- Slat Deck provides strength, no-sag mattress support, no links to connect, easy to clean
- Patient weight capacity is 350 lbs. Total weight capacity 450 lbs. evenly distributed. Includes occupant, mattress, side rails, trapeze, accessories, personal affects, etc.
- Ultra light head and foot frame weigh only 41 lbs. each making one person transport and assembly easy
- #15571 Competitor II Crank Hi/Lo Hgt Adj. Bed frame only with casters
- Full and Half Rails meet IEC Entrapment guidelines when used on the #15561/15571 beds and are backward compatible for use with the #15560 & 15570 Competitor beds
- Product Safety Standards: Meets IEC 60601-2-52 60601-1 60601-1-2 60601-1-11

PRODUCTS	Description
15571	Compet II Hi/Lo Hgt Adj Bed, 1/cs
15571FR	Compet II Hi/Lo HgtAdjBed w/FRail, 1/cs
15571FR-PKG	Compet II HiLoHgtAdjBedw/FRail,Matt,1/c
15571FR-PKG-FM	CompetII HiLoHgtAdjBedw/FRail,FMatt,1/c
15571FR-PKG-TM	CompetII HiLoHgtAdjBedw/FRail,TMatt,1/c
15571HR	Compet II Hi/LoHgt Adj Bed w/HRail, 1/c
15571HR-PKG	CompetII Hi/LoHgtAdjBedw/HRail,Matt,1/c
15571HR-PKG-FM	CompetII HiLoHgtAdjBedw/HRail,FMatt,1/c
15571HR-PKG-TM	CompetII HiLoHgtAdjBedw/HRail,TMatt,1/c

SPECIFICATIONS	Description
Weight Capacity	350 lbs (158.8 kg) Patient Wgt 450 lbs (204 kg) Total Wgt Includes occupant
Carton Shipping Weight	163 lbs. (74 kg)
Weight	146 lbs. (66.2 kg)

ACCESSORIES	Description
15571TRANS	Comp II Transporter Only 1/EA, EA
15561EXT	Compet II 4" Extension 1/EA, EA

ACCESSORIES**Description**

15561TRANSC	Comp II Transporter w/4 casters 1/CV, CS
13009TRAPBV	Trapeze, Bar, Clamp, Brown Vein, 1/cv, CS